


Ancillary ECMO Products

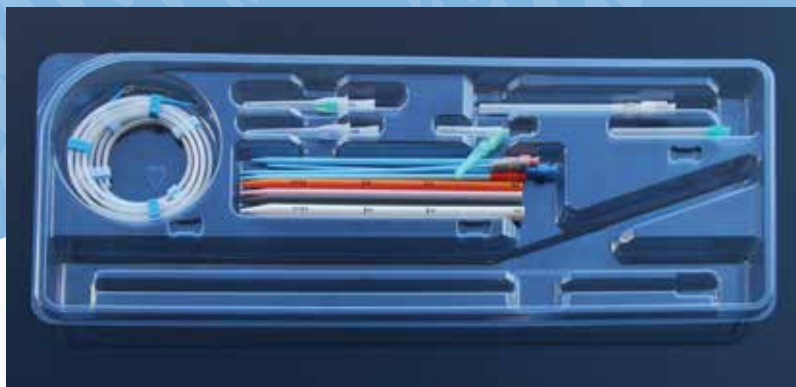

Percutaneous Introducer Kit

Part: PCT

The PCT is designed to place the OriGen ECMO catheter percutaneously to minimize the cut-down on the jugular vein during cannulation. The needles, guidewires, and dilators are all matched to the size catheter that you are using during the procedure.

Percutaneous cannulation can reduce the need for post-ECMO wound treatment and often the site wound can be treated with a simple bandage instead of sutures.

Introducer sizes: 8F, 11F, 14F, 16F, 19F, 23F, 28F, 32F

PCT for use with VV13F, VV16F, VV19F, VV23F, VV28F, VV32F

Americas & Asia

7000 Burlison Rd, Bldg D
Austin, TX USA 78744-3202
Tel: +1 (512) 474 7278
Tel: Toll Free (800) 233 9014
Fax: Toll Free (888) 812 8411 or
+1 (512) 617 1503

Europe & Middle East

OriGen Biomedical GmbH
Romerstrasse 14
72393 Burladingen, Germany
Tel: +49-7475-915591
Fax: +49-7475-914718
E-mail: sales.europe@origen.com



For availability in your country, please contact your local distributor.

Hemofiltration Tubing Set

Part: HFT

The transfer set is intended for draining fluids from the body. A set of three color coded tubes are to be connected to the correct line during ECMO. The tube set includes a venous drainage line (blue), an arterial drainage line (red), and an additional drain line (clear). Both arterial and venous lines have a sampling line with a side clamp. The drain line has a standard spike port if aseptic connection is desired.

OriGen
BIOMEDICAL

www.origen.com