

CRYOSTORE 25 FREEZING BAG

The OriGen Biomedical® CryoStore 25 Freezing Bag is intended for freezing blood components.



It provides you with the traceability and durability that you need.

ALIQUOTS - MATCHED BAG AND TUBE PRINTING

The CryoStore 25 Freezing Bag has a printed number on each chamber for easy identification. The numbers on the chambers match the pre-marked printed segments, which hold approximately 100µl of cell product. This can aid in tracing aliquot samples to the bag.

DURABILITY- BETTER MATERIAL, PROVEN PERFORMANCE

The OriGen CryoStore Bag is the industry-preferred cryogenic freezing bag. The hermetically sealed spike ports and film used in every CryoStore bag are made from a robust EVA.

NEEDLE-FREE INJECTION PORT

The “N” tubing set has luer-actuated needle-free injection ports. These swabable ports maintain sterility of the sample through numerous actuations (white paper available) and can eliminate the possibility of needle sticks.

EASE OF USE - BUILT-IN

The donor tube is intentionally placed in the corner of the bag for easier debubbling and sealing. OriGen can provide the CryoStore 25 Freezing Bag with one of several stock or special order tube sets to fit your processing needs. The tube sets may include SCD weldable tubing, needle-free access, male luer connections, and/or female luer connections.

EFFICIENCY

The chambers are filled from a single location, and the channel located at the bottom and top of the bag allows for easy heat seals between each chamber with a standard tube sealer. The separation of the individual chambers allows for easy separation of the cell product and quality sample.


OUR QUALITY STANDARDS

A Certificate of Conformance (CoC) is included with each CryoStore 25 Freezing Bag. OriGen's CryoStore 25 Freezing Bags are CE marked and FDA 510(k) cleared. OriGen is certified to ISO 13485:2016, MDSAP certified, and all products are manufactured to GMP guidelines.



CRYOSTORE 25 FREEZING BAG

The OriGen Biomedical® CryoStore 25 Freezing Bag is intended for freezing blood components.

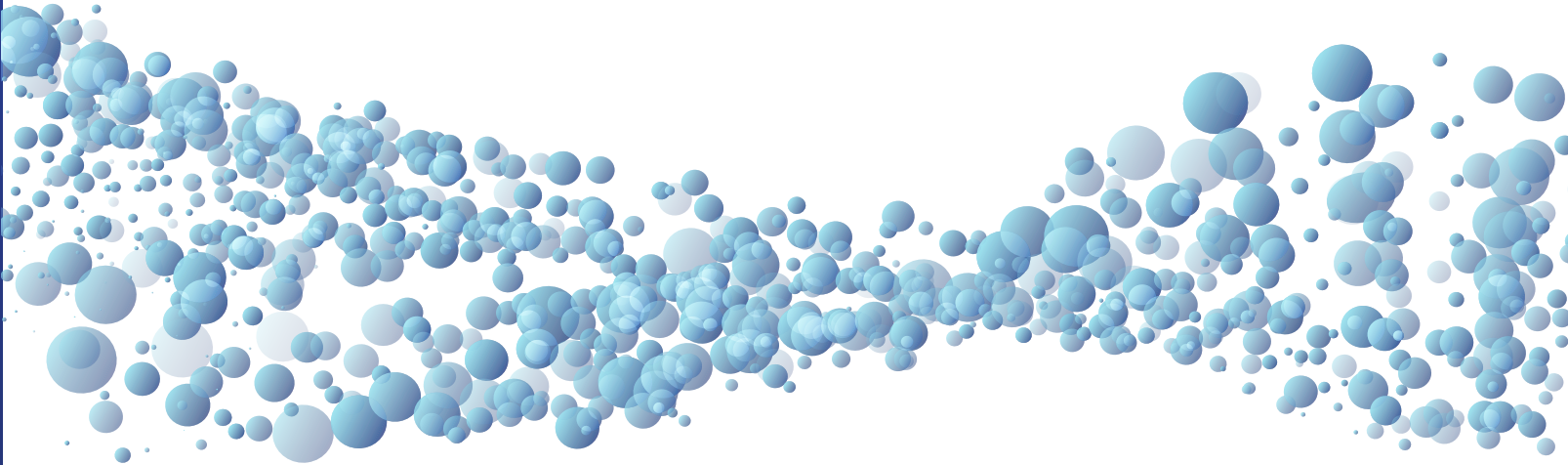

PART CODE	TUBE SET DESCRIPTION	PACKAGING	OVERWRAP	DOUBLE POUCHED	CASE QUANTITY	SPECIAL ORDER
CRYOSTORE 25 EVA FREEZING BAG, 10-25 ML						
CS25F	"F" Tube Set - Female luer with caps	Pouch of 4			24	
CS25FP	"N" Tube Set - One female luer with cap, two female luers with caps	Pouch of 1			12	✓
CS25N	"N" Tube Set - One male luer with cap, two female luers with caps, and a needle-free injection port	Pouch of 4			24	
CS25NS	"N" Tube Set - One male luer with cap, two female luers with caps, and a needle-free injection port	Pouch of 1		✓	12	
CS25NSW	"N" Tube Set - One male luer with cap, two female luers with caps, and a needle-free injection port	Pouch of 1	✓	✓	12	
CS25XN	"N" Tube Set, 8" donor tube - One male luer with cap, two female luers with caps, and a needle-free injection port	Pouch of 4			24	✓
CS25Y4	"Y4" Tube Set - One spike, one male luer with cap and 50cm SCD Tube with a needle-free injection port and pinch clamp	Pouch of 4			24	✓

7000 Burleson Rd, Bldg D | Austin, Texas 78744 | 512.474.7278 | Toll Free 800.233.9014 | sales.us@origen.com | origen.com

CE Marked | FDA Cleared | ISO 13485 Certified | MDSAP Certified | Manufactured to GMP Guidelines

ML55.R01:2023-4

© 2023 OriGen Biomedical. All rights reserved.
OriGen Biomedical® is a registered trademark of OriGen Biomedical, Inc