


# CRYOSTORE 25 FREEZING BAG

The OriGen Biomedical® CryoStore 25 Freezing Bag is intended for freezing blood components.

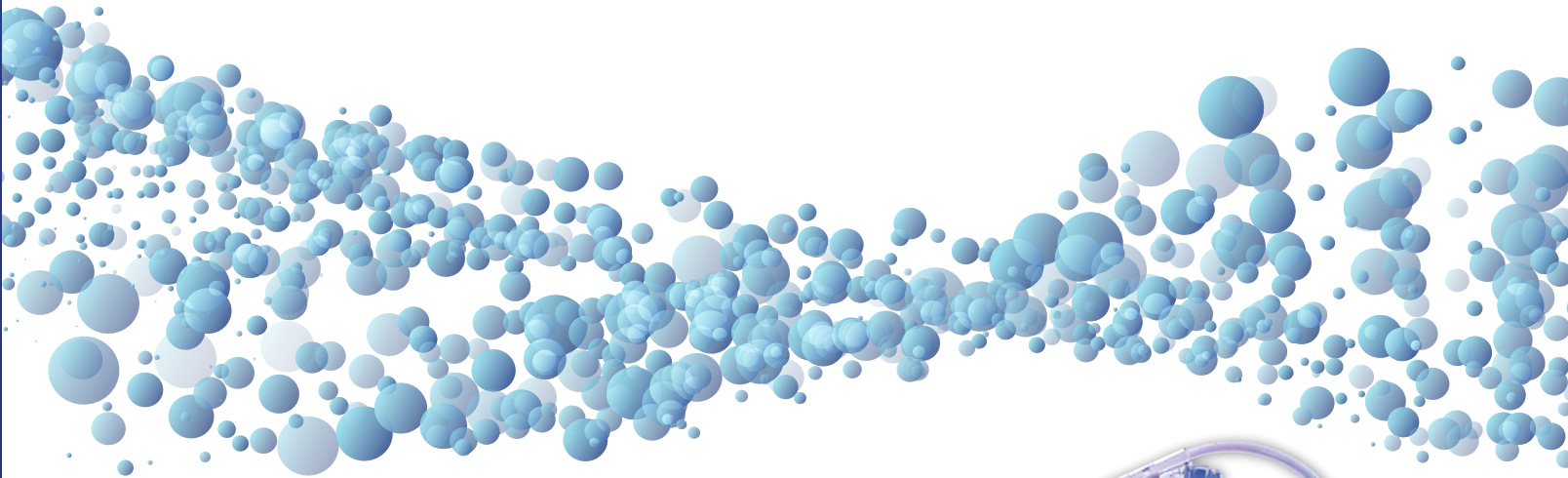

*It provides you with the traceability and durability that you need.*

## ALIQUOTS – MATCHED BAG AND TUBE PRINTING

The **CRYOSTORE** 25 Freezing Bag has a printed number on each chamber for easy identification. The numbers on the chambers match the pre-marked printed segments, which hold approximately 100µl of cell product. This can aid in tracing aliquot samples to the bag.

**DURABILITY– BETTER MATERIAL, PROVEN PERFORMANCE** The OriGen Biomedical **CRYOSTORE** Bag is the industry-preferred cryogenic freezing bag. The hermetically sealed spike ports and film used in every **CRYOSTORE** bag are made from a robust EVA.

## NEEDLE-FREE INJECTION PORT

The “N” tubing set has luer-actuated needle-free injection ports. These swabable ports maintain sterility of the sample through numerous actuations (white paper available) and can eliminate the possibility of needle sticks.

## EASE OF USE – BUILT-IN

The donor tube is intentionally placed in the corner of the bag for easier debubbling and sealing. OriGen can provide the **CRYOSTORE** 25 Freezing Bag with one of several stock or special order tube sets to fit your processing needs. The tube sets may include SCD weldable tubing, needle-free access, male luer connections, and/or female luer connections.

## EFFICIENCY

The chambers are filled from a single location, and the channel located at the bottom and top of the bag allows for easy heat seals between each chamber with a standard tube sealer. The separation of the individual chambers allows for easy separation of the cell product and quality sample.


## OUR QUALITY STANDARDS

A Certificate of Conformance (CoC) is included with each **CRYOSTORE** 25 Freezing Bag. OriGen Biomedical **CRYOSTORE** 25 Freezing Bags are CE marked and FDA 510(k) cleared. OriGen is certified to ISO 13485:2016, MDSAP certified, and all products are manufactured to GMP guidelines.



# CRYOSTORE 25 FREEZING BAG

The OriGen Biomedical® CryoStore 25 Freezing Bag is intended for freezing blood components.

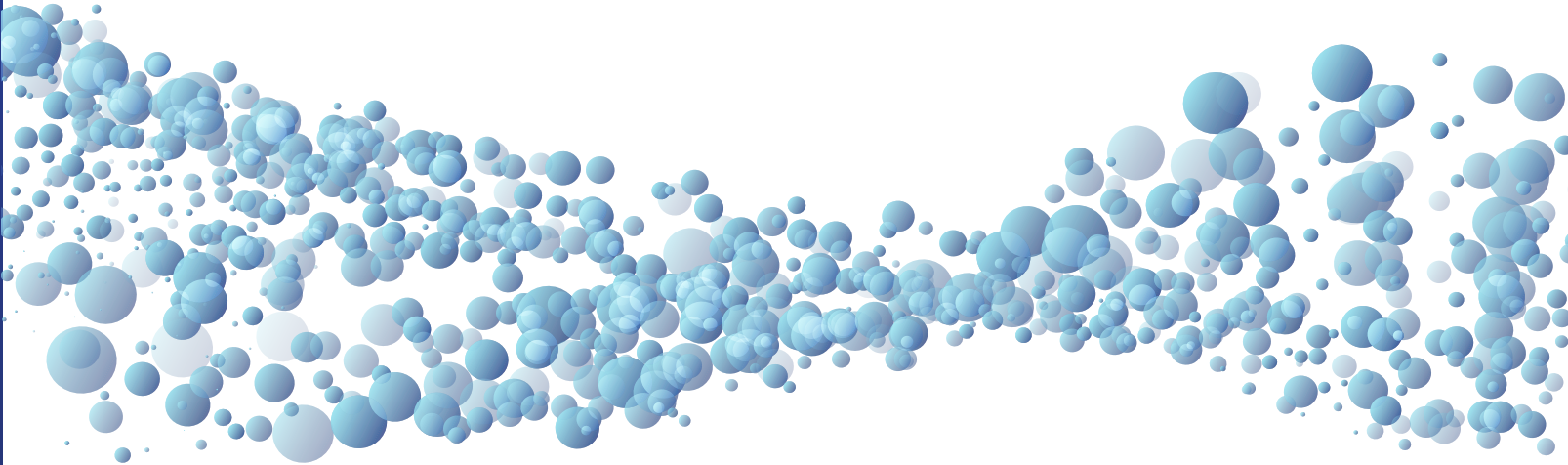

PART CODE	TUBE SET DESCRIPTION	PACKAGING	OVERWRAP	DOUBLE POUCHED	CASE QUANTITY	SPECIAL ORDER
CRYOSTORE 25 EVA FREEZING BAG, 10-25 ML						
CS25F	"F" Tube Set - Female luer with caps	Pouch of 4			24	
CS25FP	"N" Tube Set - One female luer with cap, two female luers with caps	Pouch of 1			12	✓
CS25N	"N" Tube Set - One male luer with cap, two female luers with caps, and a needle-free injection port	Pouch of 4			24	
CS25NS	"N" Tube Set - One male luer with cap, two female luers with caps, and a needle-free injection port	Pouch of 1		✓	12	
CS25NSW	"N" Tube Set - One male luer with cap, two female luers with caps, and a needle-free injection port	Pouch of 1	✓	✓	12	
CS25XN	"N" Tube Set, 8" donor tube - One male luer with cap, two female luers with caps, and a needle-free injection port	Pouch of 4			24	✓
CS25Y4	"Y4" Tube Set - One spike, one male luer with cap and 50cm SCD Tube with a needle-free injection port and pinch clamp	Pouch of 4			24	✓

7000 Burleson Rd, Bldg D | Austin, Texas 78744 | 512.474.7278 | Toll Free 800.233.9014 | sales.us@origen.com | origen.com

CE Marked | FDA Cleared | ISO 13485 Certified | MDSAP Certified | Manufactured to GMP Guidelines

ML55[01]:2023-4

© 2023 OriGen Biomedical. All rights reserved.  
OriGen Biomedical® is a registered trademark of OriGen Biomedical, Inc