


STERIZIP™ Overpouch




The OriGen **STERIZIP™** overpouch is designed to provide a clean environment for your cryogenic freezing bags when thawing in a water bath.

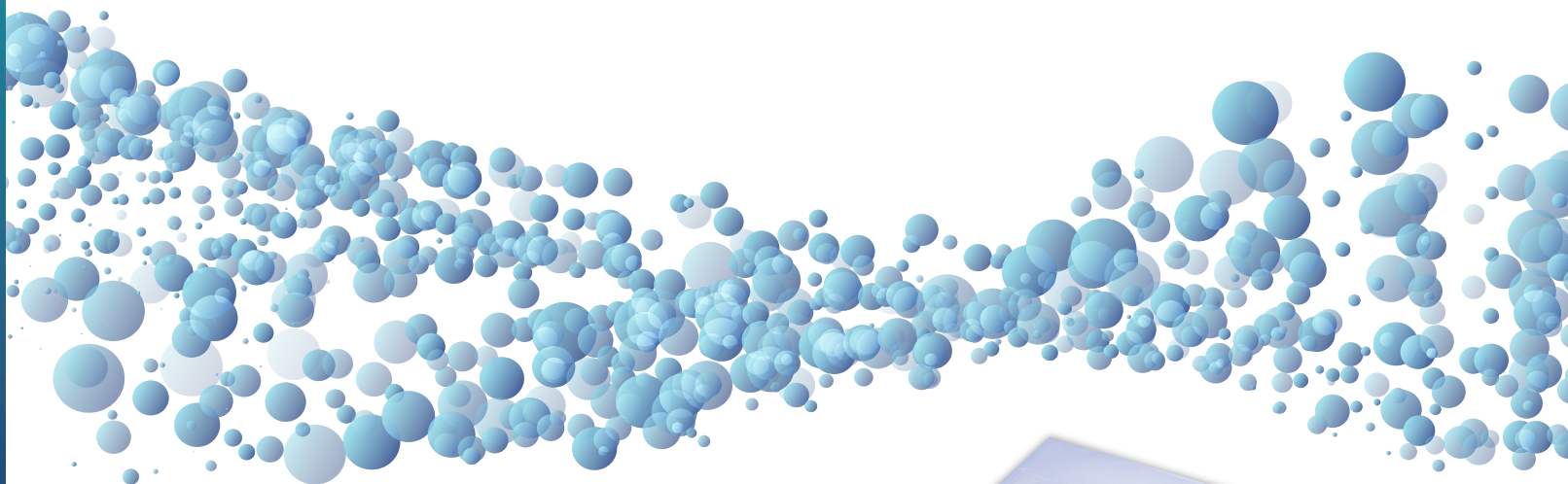

PROTECTION FOR THAWING

SteriZip™ can hold your cryogenic freezing bags during thaw and transport laboratory samples throughout your institution. Simply place your cryogenic freezing bag inside the SteriZip™ and zip the bag closed, removing air as needed. Now the SteriZip™ and contents can be placed in a water bath and massaged while thawing.

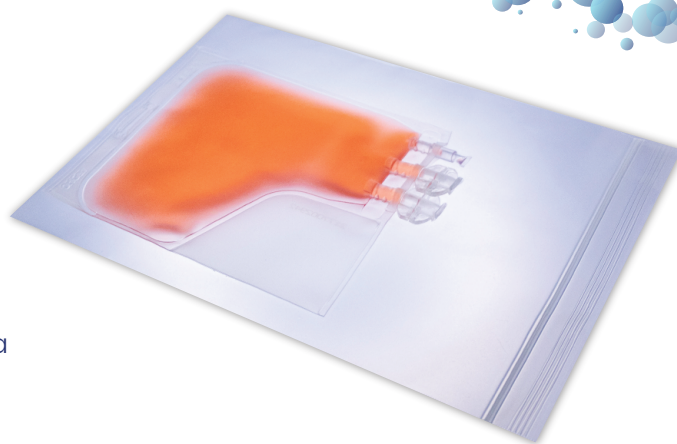

STERILE

Each SteriZip™ is inspected and pouched in our clean room, then gamma sterilized. The bags will arrive sterile.

TECHNICAL INFORMATION

PART CODE	WIDTH (cm)	LENGTH (cm)	POUCH QUANTITY	COMPATIBLE BAG SIZES
SZ1823	18	23	10	25-50ml Bags
SZ2533	25	33	10	250-750ml Bags
SZ3038	30	38	10	1000-2000ml Bags

ORDERING INFORMATION

Sterile, for single use only. To place an online order, please fill out this [form](#).

GENERAL INFORMATION

The bags are made of 0.004" thick film.